


CORONA-Fistula Probe™

Proctology

Minimal invasive laser therapy of anal fistulas

- ✓ Short fiber tip for enhanced access into smaller and curved tract
- ✓ Efficient use of laser energy for optimized closure
- ✓ Precise laser marking for guidance

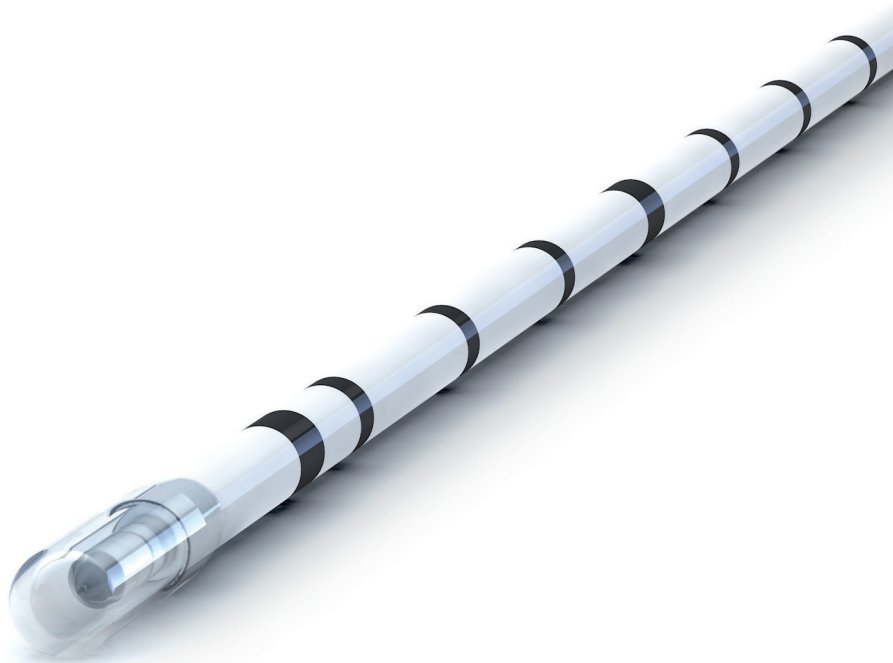

CORONA-Fistula Probe™

Proctology

The CORONA-Fistula Probe™ is specially designed to be inserted into a fistula tract to distribute the laser energy circumferentially directly to the epithelized tissue and ensures a homogenous thermal destruction of the fistula tract. The special designed short fiber tip and small outer diameter enables the reach of small and curved areas.

The dedicated laser marking allow a precise positioning of the probe within the tract.

Technical details

Outer diameter (tip)	1.0 mm
Standard length	2.5 m
Wavelength	e.g. 980 nm or 1470 nm
Typical transmission	98 %
Emission angle	60° cone angle from fiber axis
Numerical aperture	up to 0.37
Core diameter	400 µm typically

E-mail: medlaser@medlaser.srv.br
Contato: +1 866 428-2030

ORLANDO
8751 Commodity Circle Ste 5
Florida 32819

MIAMI
1001 Brickell Bay Drive
Florida 3313